

RUHOF BIOCIDE

QUATERNARY GERMICIDE

Kills All Pathogens And Viruses Including HIV-1, HBV, HCV, HAV, And Tuberculosis. Effective Against MRSA.

The safest Quaternary Germicide that is Tuberculocidal, Ruhof Biocide is basic but non-corrosive. It contains powerful quaternary germicides in a solution that is ready to use. Ruhof Biocide is a must in hospitals, clinics, nursing homes and doctors' offices. Use Ruhof Biocide to pre-clean and decontaminate critical or semi-critical medical devices prior to sterilization or high-level disinfection.

Use Ruhof Biocide to clean and disinfect hard, non-porous surfaces such as: glazed tile (bathtubs, showers, sinks), walls and furniture (examining tables, chairs, cabinets, bed springs), telephones and doorknobs. It is effective on counters and sinks and it is also a strong deodorizer for garbage cans and garbage store areas.

When used as directed, Ruhof Biocide is effective as a:

Bactericidal against Mycobacterium tuberculosis (BCG) at 20 degrees centigrade.

Pseudomonas aeruginosa, Staphylococcus aureus, Salmonella Choleraesuis, Escherichia coli, Escherichia coli O157:H7, Listeria monocytogenes, Yersinia enterocolitica, Enterococcus faecium, Corynebacterium ammoniagenes, Salmonella typhi. Ruhof Biocide meets the requirements for hospital use.

Virucidal. Ruhof Biocide, when used on hard, non-porous inanimate surfaces kills the HIV-1 Virus (Associated with AIDS), Hepatitis B Virus (HBV), Hepatitis C Virus (HCV), Hepatitis A Virus (HAV), Poliovirus Type 1, Canine Parvovirus, and Bovine Viral Diarrhoea Virus (BVDV).

Antibiotic- Resistant Strains of – Methicillin resistant Staphylococcus aureus (MRSA), Vancomycin resistant Enterococcus faecalis (VRE), Vancomycin intermediate resistant Staphylococcus aureus (VISA)


Fungicidal. Ruhof Biocide is effective against Tricophyton mentagrophytes (Athletes Foot Fungus).

Mildewstat. Controls mould and mildew on pre-cleaned, hard, non-porous surfaces.

Simply spray area to be treated with Ruhof Biocide, making sure to wet completely. Let air-dry. Repeat application at weekly intervals or when mildew growth appears.

Available in:

Cat No 345RB 12 x 500ml

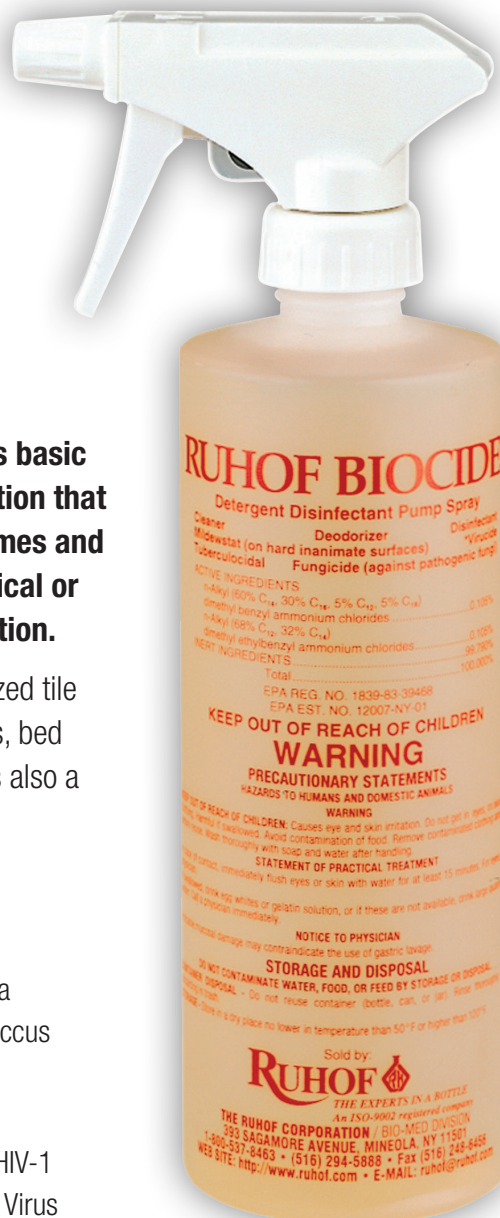

DIRECTIONS FOR USE

- Remove gross filth or heavy soil before treatment with Ruhof Biocide.
- Hold container six to eight inches from surface to be treated.
- Spray area until it is covered with the solution.
- Allow Ruhof Biocide to penetrate and remain wet for ten minutes.
- No scrubbing is necessary, wipe off with a clean cloth, mop or sponge.
- Ruhof Biocide will not leave grit or soap scum.

R **RUHOF** 
THE EXPERTS IN A BOTTLE

Peskett Solutions Ltd
Unit 4-5 Brampton Business Park
55 Brampton Road, Eastbourne, East Sussex, BN22 9AF
Tel: 01323 511038 Fax: 01323 503247
E Mail: info@ruhof.co.uk Web: www.ruhof.co.uk