



Sepsis Trolleys

Sepsis Trolleys 350

Options 351



Sepsis Trolleys

Features

- Trolleys designed for use when treating Sepsis
- Supports delivery of 'The Sepsis Six' with clearly defined locations for equipment used in each of the stages
- Door/s fitted with either:
 - Electronic (battery operated) Push Button Lock - Customer can set their own separate user, master, submaster & technician access codes
 - High Security Bolt Lock - Supplied with two keys
- Trolleys supplied with door shelves (useful storage for antibiotics etc) & an A4 document holder
 - Single Door - 2 x door shelves / 1 x A4 document holder
 - Double Door - 5 x door shelves / 1 x A4 document holder
- Six labelled drawers (Oxygen / Blood Cultures / Intravenous Antibiotics / Intravenous Fluid / Haemoglobin & Lactates / Urine Measurement)
- Optional 'Sharps Box Holder' or 'Universal Sharps Box Bracket Mounting Plate' (fits to the l/h side of the trolley)
- Sharps Box Holder - Accepts up to seven litre capacity sharps boxes (customer to fit - fixings supplied)
- Universal Sharps Box Bracket Mounting Plate - Compatible with (customer to supply the plastic brackets & fit - fixings supplied)
 - Daniels Healthcare
 - POUDS® Wall-Bracket (up to 5L) - SHARPSGUARD® mini, 1, 2.5, 3.75, 5 & com plus
 - POUDS® Wall-Bracket (7L & above) - SHARPSGUARD® 7, 11.5 & 22
 - Frontier (Sharpsafe) - 2, 3, 4 & 7
 - Sharpak Healthcare - Sharpak 18, 35, 60, 120 & 220 PLUS
- Upper work surface
- Push handles incorporate integral buffers
- 125mm polymer swivel castors
- When not in use, trolleys storing medicines should be secured to the wall. See option MT/1/6F Security Clamp

Dimensions (w x d x h)

- Single Trolley: 512 x 542 x 1072mm
- Double Trolley: 712 x 542 x 1072mm
- Drawers (internal)
 - Single Width Drawer: 338 x 297 x 115mm
 - Double Width Drawer: 338 x 297 x 115mm

Standard Colours

- Paint - Grey
- Work Surface - Orange



Electronic Push Button Lock



High Security Bolt Lock

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
SEPS/6/HSB	SEPS/6/EPB	Single Door, Six Drawer	53.4kg
SEPD/6/HSB	SEPD/6/EPB	Double Door, Six Drawer	62.2 kg

Sharps Box Holder

Features

- Accepts up to seven litre capacity sharps boxes
- Sharps box secured using a rubber band fitted with a toggle
- Attaches to the left hand side of the unit (customer to fit - fixings supplied)



Sharps box not supplied

Item	Description	Weight
SEP/SBH	Sharps Box Holder	0.38kg

Sharps Box Bracket Mounting Plate - Universal

Features

- Compatible with the following brackets (customer to supply the plastic brackets & fit - fixings supplied)
 - Daniels Healthcare
 - POUDS® Wall-Bracket (up to 5L) - SHARPSGUARD® mini, 1, 2.5, 3.75, 5 & com plus
 - POUDS® Wall-Bracket (7L & above) - SHARPSGUARD® 7, 11.5 & 22
 - Frontier (Sharpsafe) - 2, 3, 4 & 7
 - Sharpak Healthcare - Sharpak 18, 35, 60, 120 & 220 PLUS



Trolleys are compatible with most of the popular UK sharps box support brackets.

**Sharps box & bracket not supplied



Frontier (Sharpsafe)**



Sharpak Healthcare**



Daniels Healthcare Up to 5L**



Daniels Healthcare 7L & Above**

Item	Description	Weight
SEP/USBBMP	Universal Sharps Box Bracket Mounting Plate	0.16kg

Security Clamp

Features

- Allows the trolley to be secured to the 'fabric' of the building
- Requires no additional keys, utilises the trolleys own locking system
- One end is secured to the wall, the 'hook shaped' end is trapped by the door of the trolley when closed and locked
- Recommended wall fitting height 1000mm



Important Note - Security Clamps

It is a UK requirement that when the trolley is not in use it should be secured to the 'fabric of the building'. Not only does this prevent casual theft, i.e. trolleys being pushed to an isolated area and broken into, but also encourages good nursing practice, giving the trolley a 'permanent home'.


'Safe & Secure Handling of Medicines - Royal Pharmaceutical Society, GB- March 2005'

Item	Description	Dimensions	Weight
MT/1/6F	Security Clamp	Length - 950mm	0.75kg

Drawer Inserts

Features

- Plastic drawer inserts, height 100mm
- Assists segregation of drawer contents
- Suitable for single & double door trolleys
- Insert supplied with removable & repositionable dividers
 - SEP/DIS - 2 x dividers providing 4 separate areas
 - SEP/DID - 4 x dividers providing 6 separate areas
- Supplied for self-assembly
- Clip on label holders & plastic labels available to assist identification

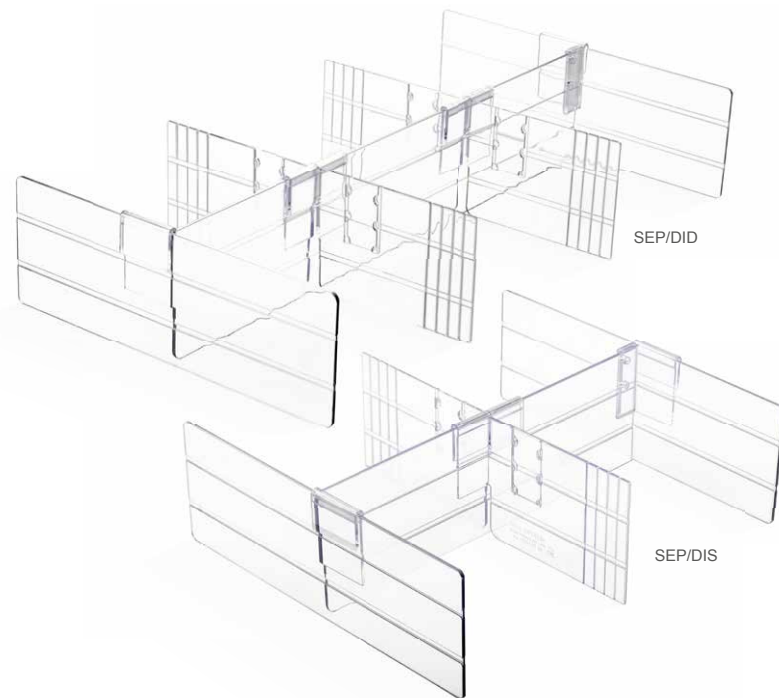


Item	Description	Weight
SEP/DIS	Drawer Insert - Set - Single Door Trolleys	0.49kg
SEP/DID	Drawer Insert - Set - Double Door Trolleys	0.55kg

Label Holders & Labels

Features

- Plastic label holders & labels
- Assists identification of drawer contents
- Label holders snap onto drawer inserts/dividers
- Accepts ISO labels 40 x 60mm (w x h)
- Sold in packs of ten

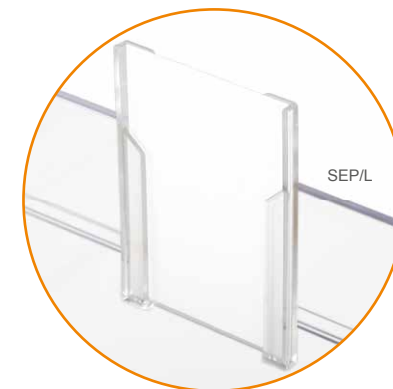

Item	Description	Weight
SEP/LH	Label Holder	0.006kg
SEP/L	Label	0.001kg

Label - Sepsis Tiger


Features

- Sammy the sepsis tiger label developed for use on Paediatric trolleys
- Compatible with both single & double door trolleys



SEP/D/6/EPB/SEP/TIG01
Trolley complete with label

Item	Description
SEP/TIG01	Label - Sammy Sepsis Tiger

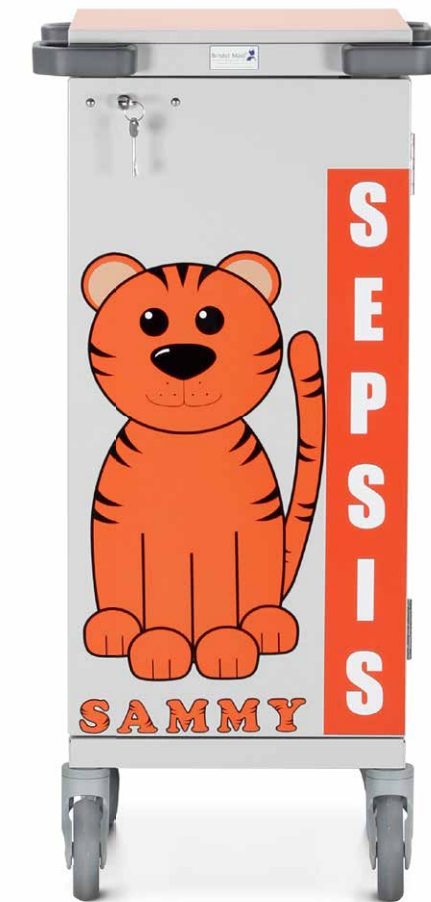

SEPS/6/HSB/SEP/TIG01
Trolley complete with label